# Lipase VET control kit

**For veterinary use only!**

For quality control procedures of the quantitative determination with the Lipase VET test kit on a solo or CUBE-VET analyser

Eurolyser Diagnostica GmbH  
Bindergasse 3  
5020 Salzburg/Austria

Tel. +43 (0)662 / 43 21 00  
Fax +43 (0)662 / 43 21 00-50  
www.eurolyser.com

## English

	Order information	Indication	Kit size
	Order number: VT 0180 Order number: VT 0181 Order number: VT 1800	Lipase VET test kit Lipase VET test kit Lipase VET control kit	16 tests/pack 6 tests/pack 1 x 5 ml (decision level)

**Test kit preparation: Allow single test at least 10 minutes to warm up to room temperature (20 - 25 °C) by placing the test into the test kit rack. Put test kit package back into refrigerator.**

### Summary

This Lipase VET control kit is a valuable tool for assuring quality and precision of the Lipase determination with the Lipase VET test kit (VT0180 / VT0181) and provides a reliable marker for monitoring the accuracy of results.

### Required Material

(not supplied with the kit)

- Pipettes
- Distilled water

### Preparation and Handling

Control material is lyophilised and must be reconstituted before use!

1. Open the vial very carefully, avoiding any loss of the lyophilized material.
2. Add exactly 5 ml of distilled water.
3. Close the vial carefully and dissolve the components completely by occasional gently swirling. Avoid the formation of foam.
4. Allow the contents to reconstitute over a period of approx. 30 minutes.
5. Analyse the control in the same way as a patient sample.

**ATTENTION:** Inaccurate reconstitution of the control and error in assay technique can cause erroneous results.

### Stability and Storage

Unopened vials are stable at 2 - 8 °C until the expiration date stated on the label.

*Stability after reconstitution:*

Stored at 20 - 25 °C:	8 hours
Stored at 2 - 8 °C:	7 days
Stored at -20 °C:	30 days

Close immediately after use. Store control tightly capped when not in use. Warning: FREEZE ONLY ONCE!

### Waste Management

Please refer to local legal requirements

### Warnings and Precautions

This control kit is for veterinary diagnostic use only!  
**DO NOT INGEST!** Exercise the normal precautions required for handling all laboratory reagents.  
As no testing method can rule out the potential risk of infection with absolute certainty, the material should be treated just as carefully as a patient specimen. In the event of exposure the directives of the responsible health authorities should be followed.

Set the species to „Control“ before running the QC material

Compare the results with the values indicated in the lot-specific value table

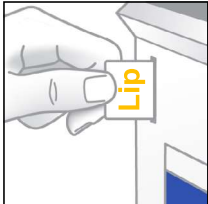
If the result is not within the acceptable range, repeat measurement. If result is again not within the acceptable range, please contact your distributor

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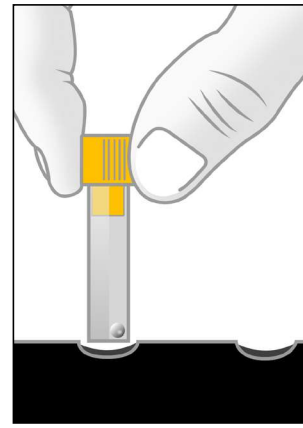
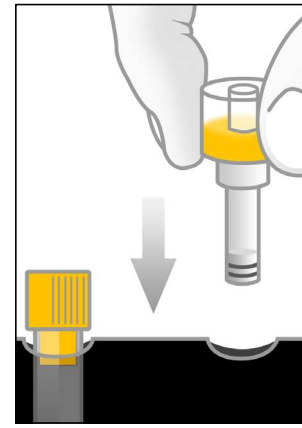
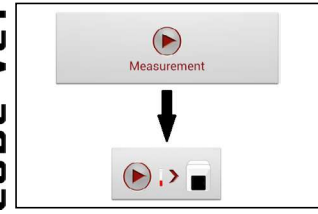
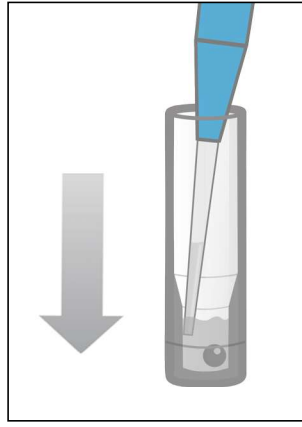
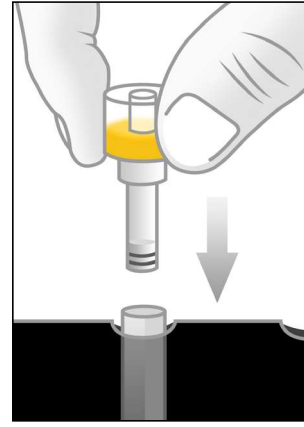
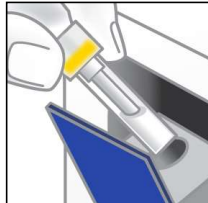
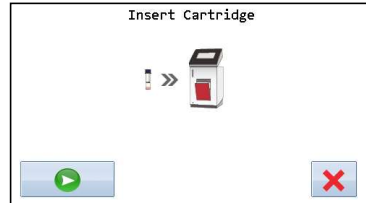
## English

### ATTENTION!

Allow single test at least 10 minutes to warm up to room temperature (20 - 25 °C) before use!

- 1. Preparation of test system**
  - 1.1 Place RFID card
  - 1.2 Place ERS cuvette in test kit rack
  - 1.3 Place ERS cap in test kit rack
  - 1.4 Press „Measurement“ button, enter required information using the touchscreen

- 2. Control preparation**

*For further details please see user manual of laboratory photometer*

  - 2.1 Set species to „Control“
  - 2.2 Set sample type to „Pipette20µ“
  - 2.3 Aspirate 20 µl control liquid out of control vial using a pipette

- 3. Control processing**
  - 3.1 Dispense 20 µl control liquid INTO THE LIQUID in the ERS cuvette
  - 3.2 Apply ERS cap firmly onto ERS cuvette
  - 3.3 Place ERS cartridge into laboratory photometer
  - 3.4 Start automatic control processing by pressing the  start button on the solo laboratory photometer, or by closing the door of the CUBE-VET laboratory photometer.

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Technical details subject to change without notice.

Doc-ID: bee4c3cb-344-f499d-87d0-9c03c6383f0f ; latest revision: DP 3.0 ; date of print: 2018-06-13

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